
SPATIAL

cod. 822333 rev. 10/2026

Air Liquide HEALTHCARE

Air Liquide Medical Systems S.r.l.
Via Bisceglie, 66 - 20152 Milano - Italy
Plant: Via dei Prati, 62 - 25073
Bovezzo (BS) ITALY
Tel. (+39) 030.201.59.11
Fax (+39) 030.2098329
medicaldevice.airliquide.com
it.medicaldevice.airliquide.com

CE

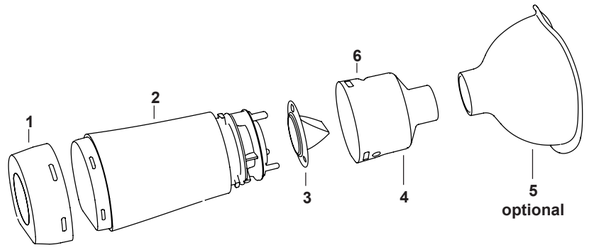ITALIANO

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ENGLISH

⚠ Read these instructions carefully before use and keep them on hand.

Intended use: This medical device belongs to the category "Valved Holding Chamber", that is intended to be used with metered dose inhalers (MDI). The device allows to improve the efficacy of drug administration at lung level by reducing side effects.

Indication for use and intended users: The device is intended for use by patients under the supervision of a doctor. Spatial is a valved holding chamber suitable for use with metered-dose inhalers (MDI) and can be used by both children and adults, at home or in hospital.

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ENGLISH

Due to the inhalation and exhalation valve, the innovative design and the antistatic material, the device improves the therapy performed with MDI.

The device is composed of:
a bottom (1) with universal opening for the MDI inhaler, holding chamber (2) inhalation and exhalation valve (3) upper body with mouthpiece (4) and air outlet (6)

The following masks are available as an optional, which can be connected directly to the mouthpiece:
- PVC pediatric mask for children between 0 and 6 years (not sterilizable. dead space: 48 ml)
- PVC adult mask for patients over 6 years old (not sterilizable. dead space: 91 ml)
- silicone infant mask for children between 0 and 2 years (sterilizable. dead space: 50 ml)
- silicone child mask for children between 2 and 6 years (sterilizable. dead space: 74 ml)
- silicone adult mask for patients over 6 years (sterilizable. dead space: 112 ml)

Technical characteristics:
Chamber dimension: 160 ml

⚠ Warning:
- Always follow your doctor's instructions regarding the type of medicine to be used and the recommended dosage.
- If you have any doubts concerning the use of this product, please consult your doctor or pharmacist.
- Do not use the product if any part is missing or damaged
- Some parts are very small and could be swallowed by children. Do not leave such parts unattended. An adult should be present at all times.
- The product must be cleaned as indicated in the following instructions before you use it for the first time, and then at least once a week.
- The repeated disinfection cycles may alter the colour or the appearance of the device. Always check the intact state of the device before use.
- It is advisable to replace the device after a year of continuous use.
- Read the instructions of the drug before use.
- There are no contraindications to the use of the device.
- Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

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Symbols

Manufacturer	Carefully read the instructions for use	Read the warnings carefully	Temperature limits	Pressure limits	Relative humidity limits				
The presence of latex was not detected	The presence of phthalates was not detected	Dispose of the device properly	CE Marking	Production date	Lot number	Medical device	Product reference	Unique Device Identifier	Swiss authorized representative

Updates to technical features
In order to continuously improve performance, safety and reliability, all medical devices produced by Air Liquide Medical Systems S.r.l. are periodically reviewed and improved. Operating instruction manuals are therefore amended to ensure constant compliance with the characteristics of the devices launched to market. If the instruction manual accompanying this device is lost, a copy corresponding to the supplied device may be obtained from the manufacturer by specifying the data shown on the technical data label.

Copyright
All information contained herein may not be used for purposes other than originally intended.
This manual is owned by Air Liquide Medical Systems S.r.l. and it may not be reproduced, either wholly or in part, without the prior written permission of Air Liquide Medical Systems S.r.l.. All rights reserved.

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ENGLISH

INSTRUCTIONS FOR USE

Before use, carefully examine the product making sure that it is assembled correctly and there are no foreign bodies or missing components. Remove any foreign objects prior to use.	Read the instructions supplied with the MDI. Remove the cover of the MDI. Insert the MDI into the bottom.	Place the mouthpiece between your lips or place the mask on your face so as to cover your nose and mouth and make sure the mask perfectly adheres to your face.	Press the MDI and inhale through your mouth, follow the breathing procedure indicated by your doctor, the exhaled air will come out from the opening (6) on the mouthpiece.

CLEANING, DISINFECTION AND STERILIZATION

Remove the bottom. Detach the mask from the mouthpiece, turn the mouthpiece and remove the valve	Cleaning: Soak and wash all parts with warm water and neutral non-abrasive cleaner. Disinfection: after cleaning, soak the parts in a disinfectant solution (e.g. a solution containing sodium hypochlorite in concentrations not greater than 2%). Follow the instructions of the manufacturer of the disinfectant solution, and make sure it is compatible with the materials of the device. Rinse thoroughly with clean water.	NB: The following procedures are not compatible with the PVC mask. The device can be washed in the dishwasher with a special detergent at temperatures not greater than 70°C, placing it at the top of the dishwasher. The parts of the device can disinfected boiling them in water for up to 10 minutes: the parts must not be in direct contact with the bottom of the pot. The parts of the device can disinfected using a steam sterilizer. The device can be sterilized in an autoclave, with cycle at 121 ° C, for a maximum of 30 minutes, for a maximum of 20 cycles.	Allow to air dry. Do not use dry cloths to dry the parts of the device.	Make sure the pieces are dry before you reassemble the device. Hook the bottom to the device's body. Turn the valve on the chamber by threading the 3 pins into the corresponding holes.	Insert the mouthpiece onto the body and turn it to align the flat side. Assemble the mask on the mouthpiece, if any. The upper part of the mask (that corresponds to the nose) should be aligned with the wording.

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ENGLISH

Materials and disposal
The device and/or its parts should be disposed of in compliance with current regulations and environmental protection laws.
Where there are no legal requirements, we recommend recycling. The table lists the materials of the various components of the device.
All materials used are phthalate and Bisphenol-A free and the presence of natural rubber has not been detected.

PART	MATERIAL
Bottom	Anti-static polypropylene + rubber
Body	Anti-static polypropylene
Valve	Silicone
Mouthpiece	Polypropylene

Conditions of use:
Temperature: from 5°C to 40°C
Relative humidity: from 15% to 93%
Atmospheric pressure: from 700 hPa to 1060 hPa

Storage conditions:
Temperature: from -25°C to 50°C
Relative humidity: from 30% to 95%
Atmospheric pressure: from 700 hPa to 1060 hPa

Warranty terms and conditions
The warranty is valid if it is accompanied by a receipt of purchase.
The manufacturer is liable for the safety, the reliability and the performance of the device if the latter is used according to the instructions for use, its intended use and if all repairs are carried out at the manufacturer or at a laboratory authorized by the same.
If the user fails to comply with the Operating Instructions and Warnings contained in this manual, the manufacturer's liability and warranty conditions shall be considered void. The warranty only covers defects caused by manufacturing, while it does not cover damages caused by normal wear and tear, any misuse or tampering.
If the defect is attributable to manufacture, the device will be repaired or replaced by the manufacturer, as it deems necessary.

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